


SurePAGE™ precast gels

High resolution & Robust separation

MES: 20+ Min @200V

MOPS: 30+ Min @200V

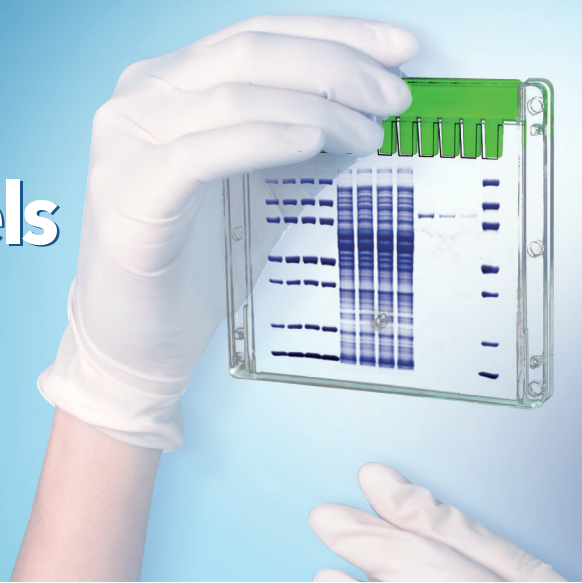

GenScript SurePAGE™ Bis-Tris gels are high-performance precast mini polyacrylamide gels that are designed to separate a wide range of protein sizes by electrophoresis. **SurePAGE gels are cast in a neutral pH buffer that minimizes polyacrylamide hydrolysis and increases gel stability. They run at neutral pH which minimizes protein modification compared to Tris-glycine gels.**

Proteins are separated on SurePAGE using MOPS or MES running buffer. We recommend using MOPS buffer for protein sizes that are >30kD and MES for proteins that are <30 kD.

We recommend using 4×LDS sample buffer (M00676) as the sample loading buffer.



Large well volume



High resolution and reproducibility



Long shelf life



Cost effective



Compatible cassette design

Superior resolution

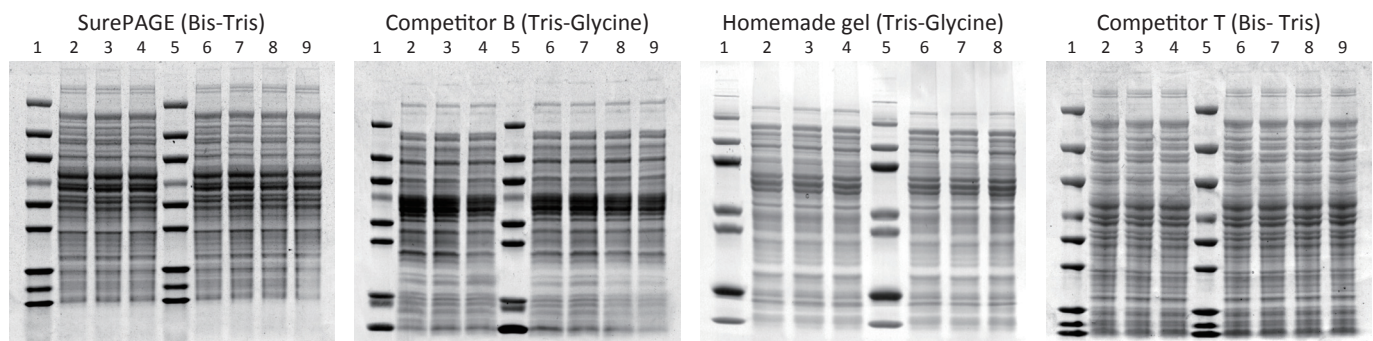


Figure 1: SurePAGE gels offer superior band resolution compared to competitors and the homemade Tris-Glycine gels.


Lane 1 and 5: protein marker (MM1397), 5 μ l. Lane 2,3,6,7,8 and 9: *E. coli* cell lysate.

Free samples (includes MOPS or MES buffer) available.
Contact us at freesample@genscript.com




Choose the gel concentration based on the migration chart below with MOPS or MES buffer.


MOPS BUFFER

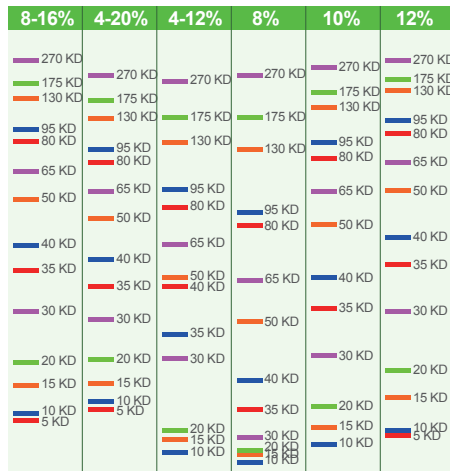

MES BUFFER

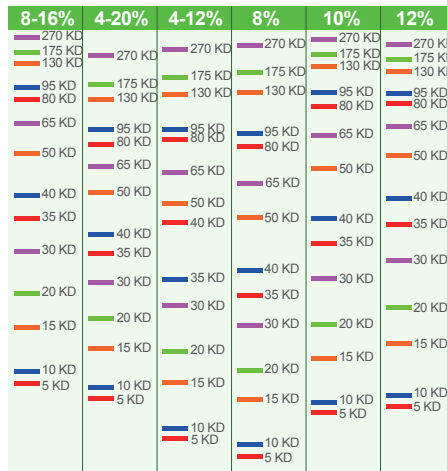

Easy to load...large well volume up to 80 **uL**

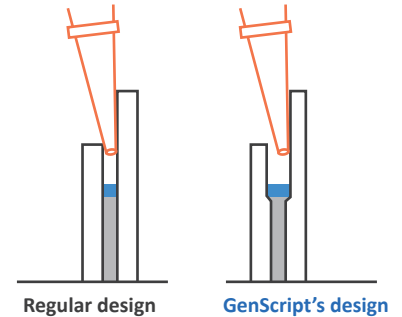


Compatible Gel Tanks

- _____
Bio-Rad Mini-PROTEAN® II & 3*
- _____
Bio-Rad Mini-PROTEAN® Tetra System*
- _____
Invitrogen Novex XCell I, II, & Surelock® (Use GenScript Gel Tank Adapter Plates)
- _____
LONZA PAGER® Minigel Chamber
- _____
Hoefler Mighty Small (SE 260/SE 250)
- _____
Hoefler Tall Mighty Small (SE 280)

*please reverse the gasket, see manual for detail

Cost effective

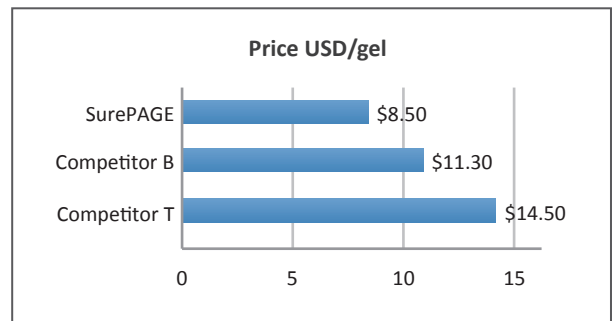

SurePAGE 10x8 gels ordering list:

% Acrylamide	Wells No.	Well volume	Quantity	Cat. No.	Price
4-12%	10, 12, 15	80, 60, 40 µl	10 gels	M00652, M00653, M00654	\$85
4-20%	10, 12, 15	80, 60, 40 µl	10 gels	M00655, M00656, M00657	\$85
8-16%	10, 12, 15	80, 60, 40 µl	10 gels	M00658, M00659, M00660	\$85
8%	10, 12, 15	80, 60, 40 µl	10 gels	M00661, M00662, M00663	\$85
10%	10, 12, 15	80, 60, 40 µl	10 gels	M00664, M00665, M00666	\$85
12%	10, 12, 15	80, 60, 40 µl	10 gels	M00667, M00668, M00669	\$85

Related products:

Product Name	Size	Cat. No.	Price
4xLDS sample buffer	10 ml, 250 ml	M00676-10, M00676-250	\$10, \$90
MES SDS Running buffer powder	5 ps/box (10 gels)	M00677	\$25
MOPS SDS Running buffer powder	5 ps/box (10 gels)	M00138	\$25
Broad Multi Color Pre-Stained Protein Standard	250 µl, 1250 µl	M00624-250, M00624-1250	\$70, \$280
GenBox Mini Electrophoresis Tank	1 unit	L00780	\$480
GenBox Mini Blot Module	1 unit	L00781	\$500