

## Bedding system for animal cages

Bedding dispenser mobile unit 

Dumping stations range 

A full range of devices for bedding management, characterized by:

- Ergonomy
- Operator safety
- Allergene containment
- Suitable for multiple bedding types



with shredder  
for vacuum system connection



mobile unit



## Steam sterilization autoclaves, LVS range

new Steelco equipments dedicated to scientific and research laboratories.

## MITO FB COMBI

Automated feeding bottles washing system




### Decapping, emptying, washing, sanitizing, filling and recapping:

*all in one solution within the smallest footprint.*

- **Over 1080 bottles processed per hour.**
- **Fully automated.** Compatible with robotic loading and unloading systems.
- **Compact solution.**  
Less than 4 m<sup>2</sup> footprint.

**Washer dimensions (WxDxH)**  
mm 3400/1150/2450 - inches 133,8/45,3/96,5

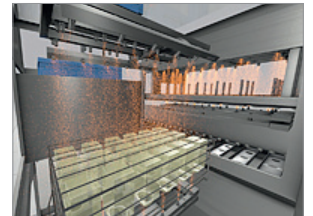

## MITO PT

automatic pass-through feeding bottle washers



- **Fully automatic process**
- **High productivity** up to 30 crates per hour.
- **Suitable for pass-through applications.**
- **Separate washing and rinsing circuits.**
- **Individual nozzles for full bottle coverage.**



Easy integration with de-capping and re-capping modules

## De-Capper module



**Modular concept units** for the automatic **de-capping** or **filling and recapping** of animal drinking bottles.

- **Easy integration** into existing washing areas.
- **Reduction of repetitive manual actions** and handling of heavy loads
- **High productivity**, up to 60 crates per hour

## Re-Capper module



## MITO FBS

semiautomatic feeding bottle washer with integrated filling system




"in line" configurations: manual de-capper, ultrasonic bath for cap treatment, dripping station

**All in one bottle treatment solution:**

- Patented rotating system for the automatic **emptying, washing, rinsing** and **filling** of bottles
- **1x18 bottle crate per cycle**
- Independent washing and rinsing system
- Full bottle coverage with individual nozzles
- Easy access to the water filtration system

**The core of the system**

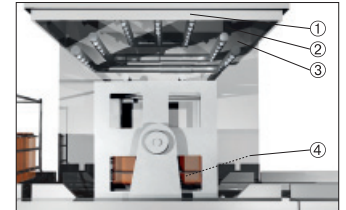

Water distribution system:

- ① filling system
- ② top rinsing system
- ③ washing system
- ④ bottom rinsing system



## MITO 2B

semiautomatic feeding bottle washer

- **2x18 bottle crates per cycle**
- **Easy integration with:** loading and unloading tables, ultrasonic cap treatment, filling modules
- Independent washing and rinsing system
- Full bottle coverage with individual nozzles



washing, rinsing and filling solution



## Water treatment solution for feeding bottles

### Filtration unit for drinking water

High level water filtering system on a wall-mounted structure.

### Acidification unit for drinking water

The system is equipped with an automatic proportional outlet pump and with two pH meters.

### MITO filling module feeding bottle filling station

Filling module for one or two bottle crates. This module can be manually loaded or automatically conveyed.

