

Wobble-not™

Low Insertion Force Serological Pipet


U.S. Pat. 9,566,579



Request a free sample pack




LOSE THE WOBBLE
LOSE THE
DRIP

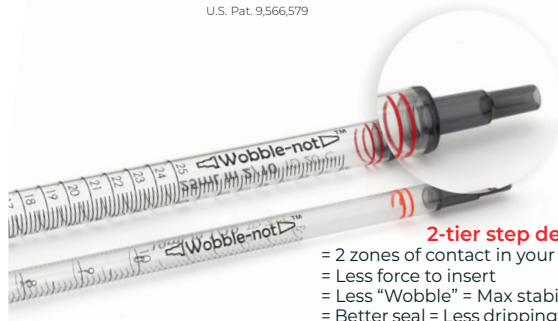

Do any of your other serological pipets look like this?

2-tier step design

- = 2 zones of contact in your pipet controller
- = Less force to insert
- = Less "Wobble" = Max stability
- = Better seal = Less dripping = Better accuracy

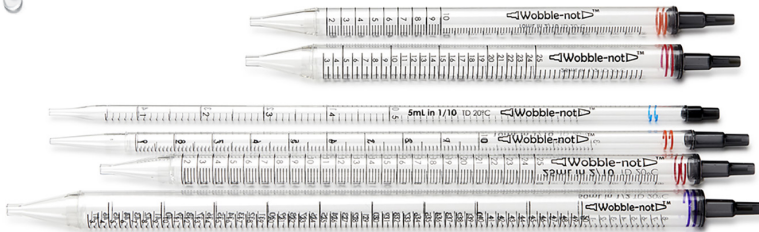

Paper/Plastic indiv. wrapped Standard Length

Volume	Product #	qty per sub-pack	qty per case
5 mL	4090-0005	50	200
10 mL	4090-0010	50	200
25 mL	4090-0025	25	200
50 mL	4090-0050	25	100

Paper/Plastic indiv. wrapped Shorty Length

Volume	Product #	qty per sub-pack	qty per case
10 mL	5090-0010	25	200
25 mL	5090-0025	25	100

Features of all Wobble-nots pipets


- Sterile - RNase/DNase free
- Non-pyrogenic
- Filter made of porous polypropylene
- Ultra-clear markings
- Accuracy $\pm 2\%$ at full volume

Bulk packaged Standard Length


Volume	Product #	qty per sub-pack	qty per case
5 mL	4290-0005	25	500
10 mL	4290-0010	25	400
25 mL	4290-0025	10	150
50 mL	4290-0050	10	100

Standard Pipet

If you don't use enough insertion force, the pipet is loosely held and is wobbly.

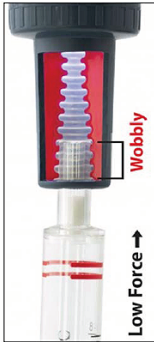

You must use higher insertion force to get less wobble and a more stable pipet.

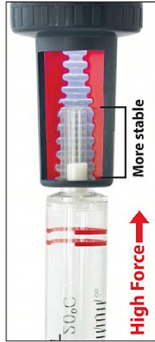

Wobble-not™


With low insertion force, a Wobble-not pipet is held very stably.



All controllers have the same design for the nozzle insert to hold onto any pipet.



Why it works

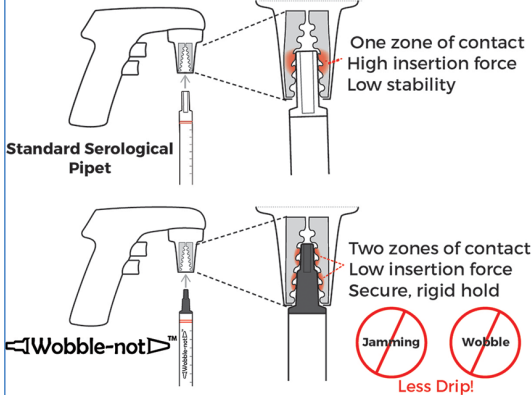

wobble; wob-ble (wōb'əl); noun.
A movement or rotation with an uneven or rocking motion or an unsteady motion from side to side.

