

Soft Dry Bath -one for all shapes of vessels

Heating in reflux, distillation and rectification, Oil free heating, completely new experience

WIGGENS launched one for all shapes of vessels Soft Dry Bath which can be easily mounted on the top of 135mm round hotplate top.

- > The high quality soft dry bath developed and produced by Wiggins have uniform heat transfer and fast heat conduction, which is very suitable for replacing the small oil bath commonly used in the laboratory dry bath module, electric heating sleeve and other heating methods.
- > Aluminum bottom of the bath, it helps to transfer the heat to alloy particles very fast. Surrounded by PTFE material, for insulation and anticorrosion
- > With infrared heating (IR) technology and ACC self-tuning temperature control technology of Wiggins, it can achieve faster heating effect and better temperature stability.


- Suitable for any shape of container
- Replacement of Silicon oil
- Replacement of thermostat module
- High temperature
- Safety
- New experience
- No cleaning required.
- Energy conservation and environmental protection
- The exclusive technology



One For All



Suitable for any shape of container

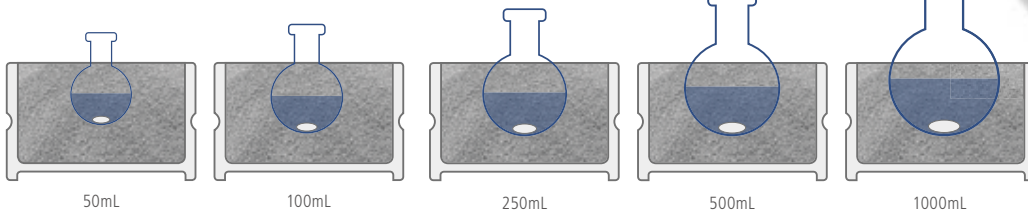

Suitable for different volume of containers

Order information

Order No.	Model	Description
SDB-1	Dry bath tank and jacket	The dry bath tank and PTFE insulation jacket internal diameter: 140mm
SDB-2	Alloy seeds	The alloy seeds, 2kg

* For digital hot plate / stirrers (WH240-R or WH260-R), you need to order separately.



For stirrer bars,
please reference
Page 73